
Atlanta★

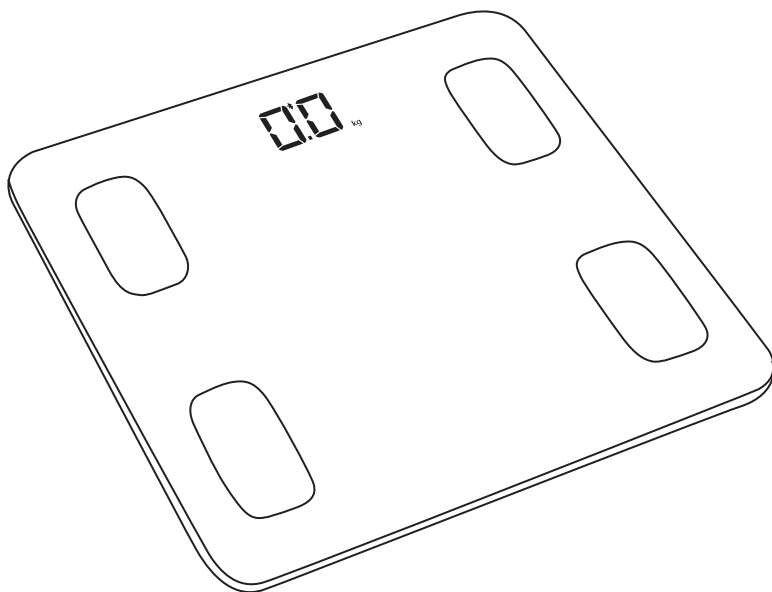ATH-6163

Personal electronic scales

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
Ver.2.06.2010

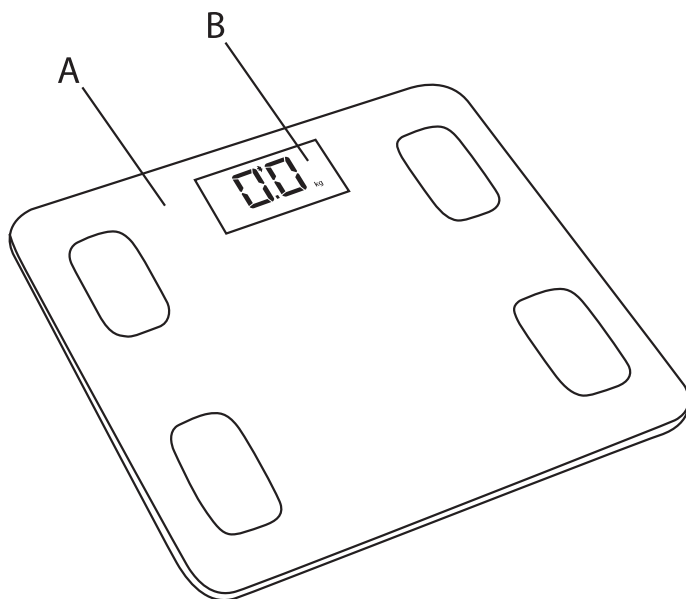
DEAR CUSTOMER!


CONGRATULATIONS ON YOUR GOOD CHOICE! NOW YOU ARE THE OWNER OF ATLANTA™ PRODUCT NOTABLE FOR ITS PROGRESSIVE DESIGN AND GOOD QUALITY PERFORMANCE. WE HOPE OUR PRODUCTS WILL STAY WITH YOU FOR YEARS ON END.

BEFORE USE

1. Please, attentively read this instruction before the beginning of operation with scales. Insert four AAA batteries as shown in the correct orientation when installed (batteries not included).

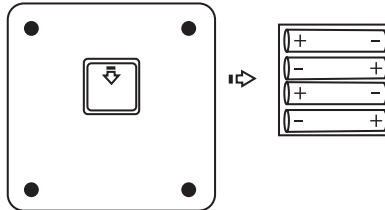2. Next to the battery compartment, there is a button for changing the weight unit. Click if you want to change the displayed unit of measurement. You will see its change on the LCD.

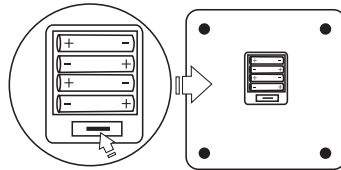3. Go to the "Google Play" or "App Store" app store on your smartphone. Install the "SenssunLife" or "MovingLife" app. You can also find apps by scanning the QR code below:

Requires iOS 8.0 or later

EN
3-9

Requires Android 4.3 or later


RUS
10-16

KZ
17-23

App: SenssunLife

App: MovingLife

Open the installed app and turn on Bluetooth on your phone. Register your account and follow the step-by-step instructions to set up your profile and connection. The scale should be turned on when you start making the connection. All measurements made can be saved in the history for detailed analysis at any time and in any place by connecting to your account.

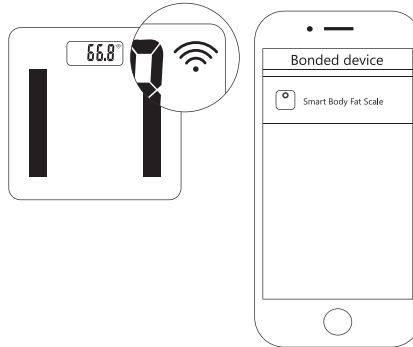4. Symbols of connection status:  or 

Fail to be connected: the symbol will be flashing

Succeed to be connected: the symbol will be on constantly

Some smartphone models may not have the characters connect, for more information you can find in the user manual of your smartphone.

SAFETY MEASURES

- Electronic scales are designed for operation from the power supply with characteristics are specified on its marking.
- THIS DEVICE MAY NOT BE USED BY PERSONS (INCLUDING CHILDREN) WITH REDUCED PHYSICAL, SENSORY OR MENTAL CAPABILITIES OR LACK OF EXPERIENCE AND KNOWLEDGE IF THEY ARE NOT SUPERVISED OR HAVE NOT BEEN INSTRUCTED CONCERNING THE USE OF THE DEVICE BY A PERSON RESPONSIBLE FOR THEIR SAFETY.
- DON'T LEAVE THE DEVICE UNATTENDED IF IT IS SWITCHED ON OR WHEN CHILDREN ARE PRESENT NEAR IT. DON'T LET THE CHILDREN PLAY WITH OR USE THE DEVICE. CHILDREN AT THE AGE OF 8 AND OVER ARE ALLOWED TO USE THE DEVICE IF THEY ARE CONSTANTLY SUPERVISED BY ADULTS.
- In order to reduce the risk of damage of the device:
 - do not open up the housing of the product;
 - do not expose the electrical part of the device to rain or moisture;
 - do not immerse the device in water or other liquid;
- The device is designed for operating indoors, for domestic use only. Commercial use of the device is prohibited.
- In any malfunctions occurs in the course of operation discontinue to use the device. Refer to the nearest service center. Never fix problems by yourself.
- Do not place the device close to any sources of heat, such as: heating radiators, gas stoves; or in areas subject to direct sunlight, excessive dust, mechanical vibration or impact loads.
- After transportation under low temperature conditions, the device can be connected to the mains only after two hours in a warm room.
- Do not use any chemically active liquids (petrol, kerosene, alcohol, various solvents) for cleaning the housing of the device.
- Any faults in storage, transportation or operation of the device are not covered by the warranty.
- It is not allowed to use the scale if there is a pacemaker.
- Pregnant women should consult a doctor before use, as there may be inaccuracies in the measurements.

EN
3-9RUS
10-16KZ
17-23

BASIC COMPONENTS OF THE SCALES

- A - Case
- B - Liquid crystal display

SPECIFICATIONS

Fat content measurement range:	5.0-50.0%
Water content measurement range:	20.0-70.0%
Analyzer measurement accuracy:	0,1%
User age:	from 10 to 85 years
User body height:	from 100 to 250 cm
Parameter memory:	12 users
Maximal weight:	180 kg
Minimum weight:	5 kg
Measurement accuracy:	50 g
Measurement interval:	50 g
Power elements:	4 x1.5V, AAA type (do not come with)
Service life:	5 years

Date of production is represented on the label of the product in the format MM.YYYY

PROPER USE OF THE SCALES

For proper operation of the analyzer, your personal data must be entered in the mobile app for scales memory, otherwise the scales will use data by default and will display correctly your weight only. It is recommended to weigh at the same time of the day, if possible (preferably in the morning), and in the same place, activity condition, then you will be able to get comparable results. Even different surface slope may dither your weight data and, as a consequence, your body content values. In addition, you should consider that weighing when in clothes, after meals or drinking may lead to improper results. Pay attention to the long-term only. Short-term deviations in weight in most cases are connected exclusively with the liquid loss. If the weight decreases only, and the share falling at the fat increases or remains stable, this means that you have lost just the liquid, for example, as a result of training session, visit of sauna or diet oriented exclusively to a fast reduction of weight. If the weight increases, and the share falling at the fat decreases or remains the same, this means by contrast that you managed to gain in precious muscular. Simultaneous reduction of the weight and the share falling at the fat means that your diet works right and your weight actually decreases.

OPERATION

1. Install the scales on a firm, enduring and horizontal surface.
2. Place the scale on a firm, stable, horizontal surface.
3. Carefully stand on the scale with your bare feet on the contact plates. You should be positioned on the scale so that you stand on the right and left electrodes with different feet. The legs should not touch anywhere in order to avoid the passage of the analyzing pulse in a shorter circle. Do not forget that the scale is equipped with an auto-off function, and if you do not stand on the scale for 10 seconds, the scale will turn off. Stand on the scales, try to stand still, your weight will be displayed in kg. For example, 90 kg 100 g


The scales are designed for weighing up to 180 kg. In case of any excess weight, the indicator will show error. You can change the weighing standard from kilograms to pounds. To do this, click on the button next to the battery compartment, or set the required units in the mobile app. If you want to see all the measurement data, open the mobile app that is paired with the scale.

4. After weighing, the scales will automatically turn off within 10 seconds.

NOTES

- The data that you will get using this device is not medical condition for reduction of your fat and water intake. If you have intention to amend significantly your food intake, please consult first your dietarian.
- The scales do not provide consistent results for persons aged under 10 or over 80, those with high temperature, edema or osteoporosis signs, with the pulse less than 60 beats per minute, using drugs for cardiovascular system treatment, being in excitatory state or persons with athletic musculation and for professional sportsmen.
- Do not store the scales in vertical position with batteries, as in such position the scales may automatically turn on and discharge the batteries.
- Do not leave batteries in the scales, if you have no intention to use them for a long time.
- If as a result of a double measurement, the scales showed different results, your true weight will be between these two values.
- In order to provide more accurate weighing, it should be carried out in the same place in the room. By doing so, try to stand straight.

EN
3-9RUS
10-16KZ
17-23

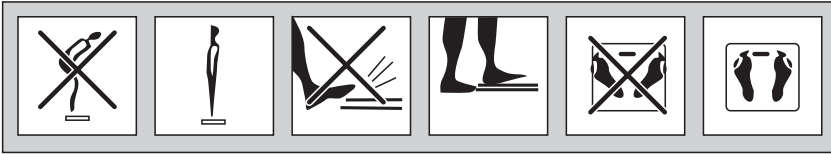- The manufacturer reserves the right to make additional alterations in the design of the electronic scales.

CLEANING AND MAINTENANCE

- The scales may be wiped using slightly moistened cloth.
- Do not use chemically active or abrasive detergents to clean the scales, as this may damage their surface.
- Do not clean the scales using tools.
- Never moisten and submerge the item in water or other liquid.

INSTRUCTIONS FOR ENVIRONMENTAL PROTECTION

The old appliances are the assemblage of technical materials and therefore cannot be utilized with household rubbish! For this reason, we would kindly ask you to support us in resource saving and environmental protection and deliver this appliance into a utilization reception center (if there is one).

TROUBLESHOOTING

Defect signs	Probable cause	Remedy
The scales do not response to button pressing or do not get activated	Electronic or static error. Low batteries	Remove batteries. Wait several seconds and then reinstall the batteries. Make sure that batteries are not low
The symbol « ____ » is not displayed	Sensor indications are out of measurement zone	Do not measure the weight of less than five kilograms
«LO»	Batteries are low	Replace batteries
Symbol «Err» or «O-Ld» is displayed	The scales were not activated by short touch-down. Maximal load of 180 kg was exceeded	Touch slightly the working surface of the scales for a short moment (activate them) and repeat the measurement. You may weigh on the scales only the weight, which does not exceed the specified maximal load. Repeat the measurement, and when measuring, stand still without motion
Symbol «Err» is displayed	Your weight is not determined, because during measurement you stood unevenly or you moved	Repeat the measurement, and when measuring, stand still without motion
Symbol «Err» is displayed	The content of the body analysed is beyond the measurement limits, your feet are in contact with each other, or poor contact of your feet with analysing plates	Repeat the measurement, having raised on the working surface of the scales with your bare feet. If necessary, moisten slightly the bottom surface of your feet. Your feet should be not in contact

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